

**ALT-ST002 – Alerta Standing Turner** - a transfer device to assist carer/attendant with transfer of the user between seat, bed, toilet, commode, wheelchair by standing and turning. Alerta Standing Turner promotes user engagement and enhances mobility.


It benefits from Alerta's GlideTurn easy swivel base, height adjustable handle, height adjustable padded knee pads, sturdy steel frame and tool free assembly.

To use this device, the user requires enough arm, leg and body strength to stand up and remain in the standing position. For users who lack these attributes, an Alerta Standard hoist or MoveAssist is recommended.

## KEY FEATURES

- > Safe Working Load: 160kg
- > GlideTurn easy swivel base
- > Height adjustable handle
- > Height adjustable padded knee pads
- > Tool free assembly and use
- > Low profile safety grip base
- > 2 x 7cm rear castors for transporting between users
- > Handle geometry facilitates nose-over-toe standing position

Do not use this device to transport a user. When patient is on board DO NOT tilt or lever the Standing Turner, the whole footplate should remain in contact with the floor at all times.

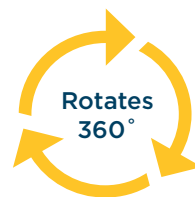

## TECHNICAL SPECIFICATION

- > Base width: 41cm
- > Handle width: 47cm
- > Handle Height Adjust: 89cm to 126cm
- > Product weight: 17.5kg